

fiocchetti
THE COLD MANUFACTURER

Fiocchetti Refrigerator



Brand Information

- Based in Italy, Made in EU
- Over 50 years history
- Dismiss the production of catering refrigerators derived
- Focus on development and manufacturing of **top quality** refrigerators for pharmacies, hospital, blood transfusion centers
- ISO 9001 and 13485 (for medical devices)
- Extend from Europe to Africa and Middle East, from South America to Australia




ECT-F Touch Features

- Touch Screen
- Triple Layer Tempered Glass
- Digital Key Lock
- Hi and Lo Temp Alarm Setting
- Door Open Alarm Setting
- Security Setting
- Built-in Temperature Logger
- Easy Recording
- SMS Function (Optional)



Alarm

- High temperature
- Low temperature
- Door open alarm
- AC power failure alarm
- Low battery for back up battery alarm

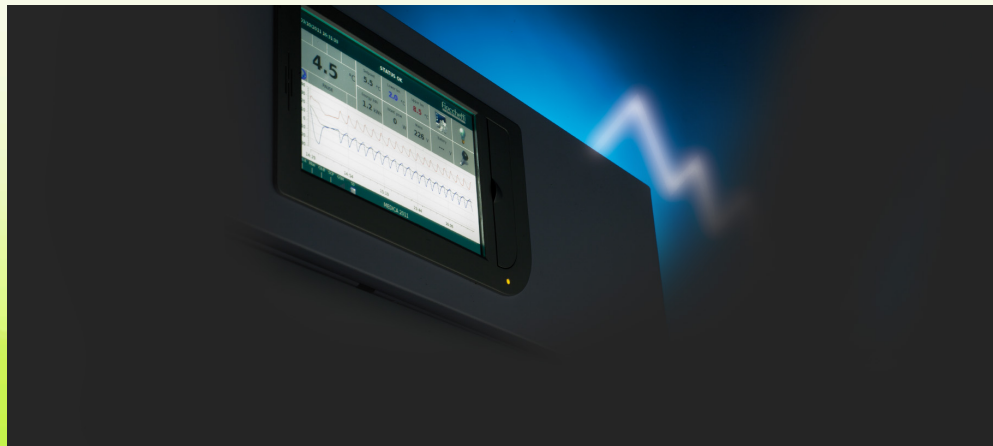

Security Setting

- Open door: by means of a gentle pressure on the key icon in the homepage
- Password protected
- Set authorize to enter different setting
- Lower chance to change critical setting
- Mechanical key supplied in the event of failure



Data Management

- Take record every 30 seconds
- Come with SD Card and can be keep record for 10 years
- USB Port for future upgrade
- SIM Card slot for SMS Function (optional)



ECT-F

- LCD Display with 4 soft touch pads
- Acoustic and visual signals for:
 - Open door
 - High Temperature
 - Low temperature
 - Power failure (when the power returns, only in event of High Temperature)
 - Dirty Condenser



ECT-F

- Record : The last 12 events are always available from the display
- Password: 2 levels of access (User and Service)
- Special functions: Stop-On-Door, Second Cooling and Advanced Defrosting allow the complete management of the refrigeration cycle, reducing temperature peaks and consumption.

Defrost

- Automatic with automatic evaporation of the condensate water

In accordance with

- 2006/42/CE, 2014/30/CE, 2014/35/CE, 2011/65/CE (RoHS) and following modifications